

Bettersize2000S喷雾激光粒度分析仪

百分页分析报告

测试范围: 0.10µm - 2000µm

样品名称:0.8mm-110kpa-2 样品编号:0001 样品来源:燕山大学 测试人员:李岳 测试时间:2024-06-016:17:25 测试单位:嘎哈医疗

物质折射率:1.333 物质吸收率:0 介质折射率:1


光学模式: Mie 分析模式: 自适应 分布类型: 体积分布

取样方法: 雾化 循环速度: 超声分散: 介质名称: 空气 分散剂: 备注:

体积平均径D[4,3] 4.886 μm 面积平均径D[3,2] 4.252 μm 数量平均径D[1,0] 3.555 μm 峰值粒径: 3.667 μm

 $D03 = 2.683 \ \mu m$ $D06 = 2.816 \ \mu m$ $D10 = 2.979 \ \mu m$ $D16 = 3.170 \ \mu m$ $D25 = 3.422 \ \mu m$

 $D75 = 5.706 \ \mu m$ $D50 = 4.186 \ \mu m$ $D84 = 6.942 \ \mu m$ $D90 = 8.059 \ \mu m$ $D97 = 10.25 \ \mu m$

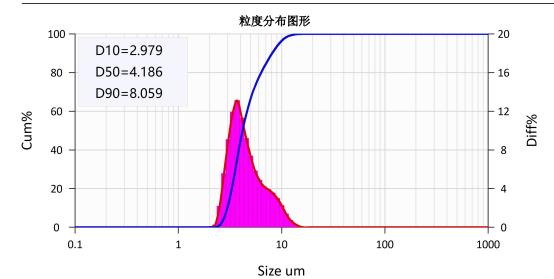

典型粒度分布表

粒径μm	含量%			
0.020	0.00			
0.050	0.00			
0.100	0.00			
0.200	0.00			
0.500	0.00			
1.000	0.00			
2.000	0.00			
5.000	66.97			
10.00	96.38			
20.00	100.00			

百分数分级表

百分位点	粒径µm	百分位点	粒径µm	百分位点	粒径µm	百分位点	粒径µm	百分位点	粒径µm
1.00	2.466	21.00	3.310	41.00	3.871	61.00	4.671	81.00	6.466
2.00	2.574	22.00	3.338	42.00	3.899	62.00	4.725	82.00	6.625
3.00	2.683	23.00	3.366	43.00	3.933	63.00	4.780	83.00	6.784
4.00	2.735	24.00	3.394	44.00	3.970	64.00	4.835	84.00	6.942
5.00	2.776	25.00	3.422	45.00	4.006	65.00	4.889	85.00	7.101
6.00	2.816	26.00	3.450	46.00	4.042	66.00	4.944	86.00	7.272
7.00	2.857	27.00	3.478	47.00	4.078	67.00	5.002	87.00	7.468
8.00	2.898	28.00	3.506	48.00	4.114	68.00	5.088	88.00	7.665
9.00	2.938	29.00	3.534	49.00	4.150	69.00	5.173	89.00	7.862
10.00	2.979	30.00	3.562	50.00	4.186	70.00	5.259	90.00	8.059
11.00	3.020	31.00	3.590	51.00	4.223	71.00	5.344	91.00	8.292
12.00	3.058	32.00	3.618	52.00	4.259	72.00	5.430	92.00	8.556
13.00	3.086	33.00	3.646	53.00	4.295	73.00	5.515	93.00	8.821
14.00	3.114	34.00	3.674	54.00	4.331	74.00	5.600	94.00	9.085
15.00	3.142	35.00	3.702	55.00	4.367	75.00	5.706	95.00	9.426
16.00	3.170	36.00	3.730	56.00	4.403	76.00	5.829	96.00	9.841
17.00	3.198	37.00	3.758	57.00	4.453	77.00	5.953	97.00	10.25
18.00	3.226	38.00	3.787	58.00	4.507	78.00	6.076	98.00	10.91
19.00	3.254	39.00	3.815	59.00	4.562	79.00	6.200	99.00	11.71
20.00	3.282	40.00	3.843	60.00	4.616	80.00	6.323	100.00	14.95

仪器制造商:丹东百特仪器有限公司 网址:www.bettersize.com 邮箱:services@bettersize.com 电话:400-655-8837